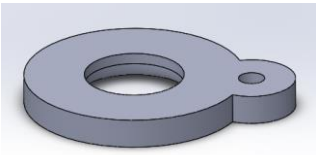
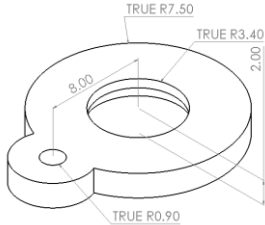
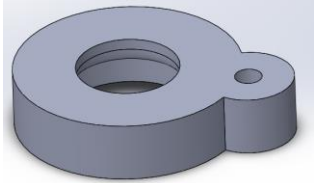
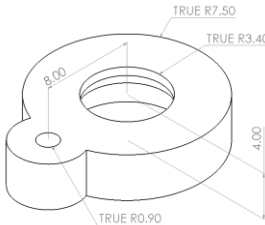
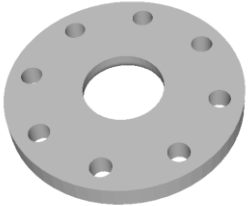
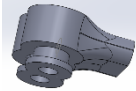
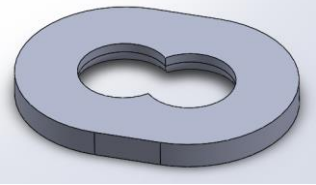
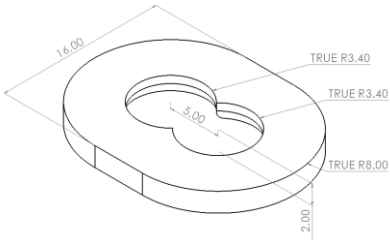
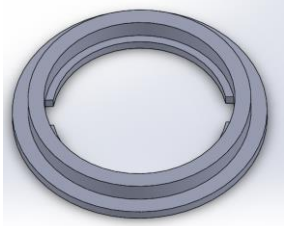
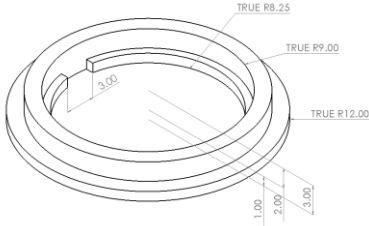
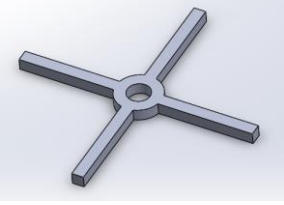
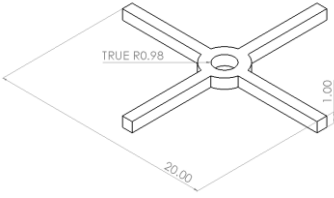

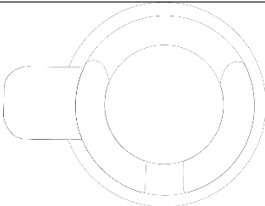
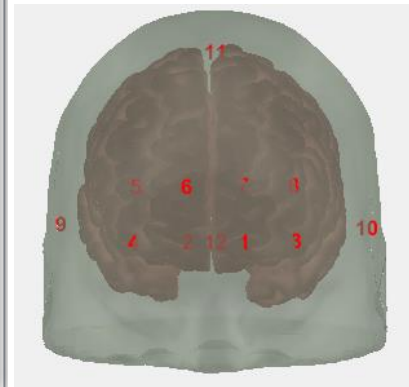
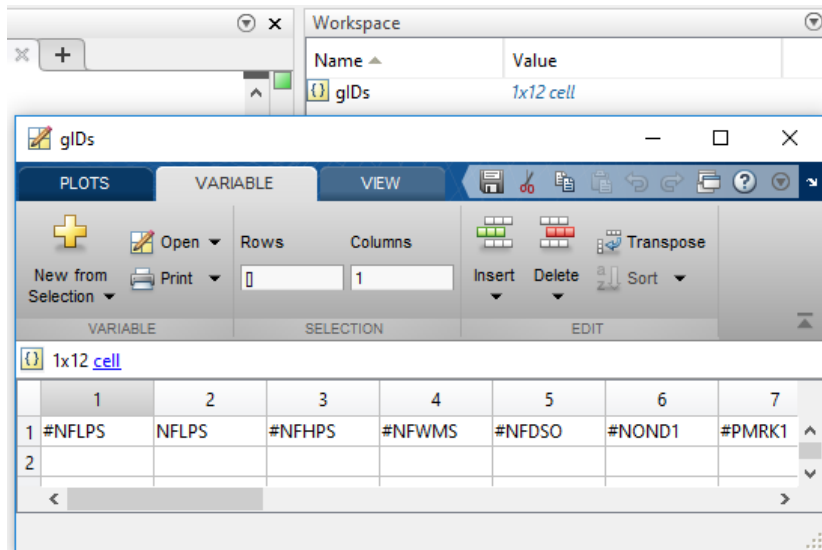


Identifier	Info	Image	Dimensions
#NFLPS	NIRS: Fiber-based optode holder (CW6 fibers) with Low Profile and single Short-separation hole		
#NFHPS	NIRS: Fiber-based optode holder (CW6 fibers) with High Profile and single Short-separation hole		
#NFWMS	NIRS: Fiber-based Wide optode holder (CW6 fibers) with Multiple Short-separation holes		Inner radius: 03.4 mm Outer radius: 08.0 mm Height: 02.0 mm SS Distance: 08.0 mm SS radius: 0.90 mm
#NFDSO	NIRS: Fiber-based Dual Ss Optode holder 		
#NOND1	NIRS: Open fNIRS Dual optode holder for optodes in printed case		
#PMRK1	Position MaRKer in small cross shape, symmetric		
#EBPAS	EEG electrode holder from Brain Products: ActiCap Snap		
#EECEH	EEG Easy Cap Electrode Holder : Ring for commercial electrode holders	Tbd	Tbd
#DUMMY	Dummy Optode, nothing placed	None	None

How to provide grommet information:

Supply matlab file “gIDs.mat” with variable “gIDs” together with the AtlasViewer.mat file in the main ninjaCap folder.

“gIDs” is 1 x N matlab cell array of strings, where N is the number of all probes in the AtlasViewer patch including dummy (anchor) optodes. The type of each element to be placed is identified by the string ID from the above table and its index. The index corresponds to the optode number in AtlasViewer.



STL files for the elements are placed under
.../ninjaCap/stl/elements/<identifier>/...

where <identifier> is the folder name identical to the defined string ID. Files within are always “grommet.stl” and “grommetComplement.stl”, where “grommet” is the element to be placed, and “grommetComplement” is the keep out geometry used to create corresponding holes/spaces in the panel grids.